

EVERYDAY

# Freedom Maverick™ Comfort AT

## Stride Comfortably!

- Low build height for those with clearance challenges
- Split keel for more ground compliance and versatility
- Increased durability without compromising the desired energy storage/return
- Smooth roll-over during stance phase of gait
- Lightweight fiberglass
- Innovative heel/keel design
- Waterproof




Split Toe



Sandal Toe

Maverick Comfort AT (F22)

# Freedom Maverick™ Comfort AT

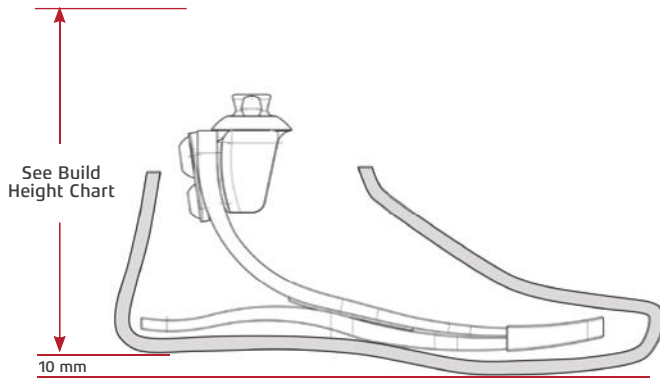

Note: Build height measurements represent a fully assembled unit.

**Product ID:** F22

**Sizes:** 23-30 cm

Sandal toe available for sizes 23-28 cm\*

**User Weight Rating:** 166 kg

**Heel Height:** 10 mm

**Connector:** Pyramid

**Warranty:** Foot module 36 months,  
Cosmetic foot shell 6 months

\*When ordering the Sandal Toe option, please contact Customer Service for lead time

## Accessories

- Foot Shell Removal Tool: (ACC-00-103)
- Spectra Sock: (ACC-00-101xx)  
30 cm, 35 cm, 40 cm
- Stiffening Bumpers:
  - 23-25 cm (KIT-00-11433)
  - 26-28 cm (KIT-00-11434)
  - 29-31 cm (KIT-00-11435)

Note: Spectra socks and heel bumpers are included in the delivery.

## Build Height & Weight

Size (cm)	Height (mm)	Weight (g)
23-25	116	626
26-28	117	763
29-30	119	841

Product weight includes foot shell. Measurements are approximate and will vary due to category selection and manufacturing variability.

## Stiffness Category Selection

Impact Level			Weight (kg)
Low	Moderate	High	
1	1	2	44-52
1	2	3	53-59
2	3	4	60-68
3	4	5	69-77
4	5	6	78-88
5	6	7	89-100
6	7	8	101-116
7	8	9	117-130
8	9		131-147
9			148-166

## Complete Prosthetic Foot (incl. foot shell) Part Number

**F22-00-1**  **A**  -  -

Category: 1-9    Size: 23-30    Version: A-mid profile, S-sandal toe    Side: L-left, R-right    Generation: 3-regular toe, 2-sandal toe    Skin Tone: L-light, M-medium, D-dark

### Foot Module Part Number

**F22-00-0**  **A**  -

Category: 1-9    Size: 23-30    Toe Configuration: RU-regular toe, Universal, SL- Sandal toe, Left, SR- Sandal toe, Right

### Foot Shell Part Number

**FTC-**  **A-1**   **4-**

Generation: 2-sandal toe, 3-mid profile    Size: 23-30    Skin Tone: L-light, M-medium, D-dark    Toe: R-regular, S-sandal    Side: L-left, R-right

See the Cosmetic Foot Shell page for more detailed ordering instructions.