

EVERYDAY

Freedom Maverick™ Xtreme AT

**Xtreme energy return, Xtreme durability,
Xtreme flexibility**

- Split keel/heel design provides excellent inversion and eversion for user stability on uneven terrains
- High impact, extreme user
- Waterproof
- Increased durability
- Natural and flexible feeling
- Lightweight fiberglass option




Split Toe



Sandal Toe

Maverick Xtreme AT (F21)

Freedom Maverick™ Xtreme AT

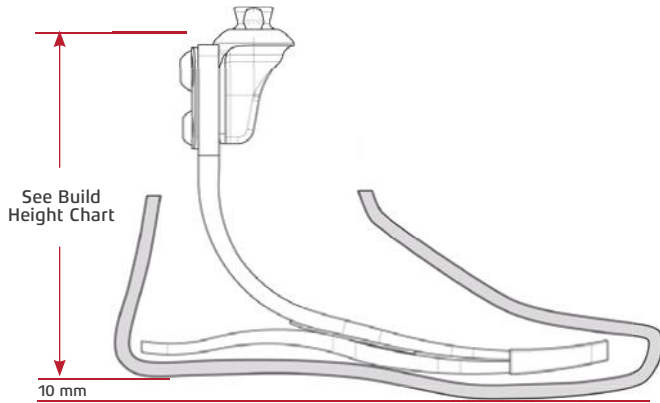

Note: Build height measurements represent a fully assembled unit.

Product ID: F21

Sizes: 23-30 cm

Sandal toe available for sizes 23-28 cm*

User Weight Rating: 166 kg

Heel Height: 10 mm

Connector: Pyramid

Warranty: Foot module 36 months,
Cosmetic foot shell 6 months

*When ordering the Sandal Toe option, please contact Customer Service for lead time

Build Height & Weight

Size (cm)	Height (mm)	Weight (g)
23-25	162	646
26-28	170	798
29-30	183	891

Product weight includes foot shell. Measurements are approximate and will vary due to category selection and manufacturing variability.

Stiffness Category Selection

Impact Level			Weight (kg)
Low	Moderate	High	
1	1	2	44-52
1	2	3	53-59
2	3	4	60-68
3	4	5	69-77
4	5	6	78-88
5	6	7	89-100
6	7	8	101-116
7	8	9	117-130
8	9		131-147
9			148-166

Accessories

- Foot Shell Removal Tool: (ACC-00-103)
- Spectra Sock: (ACC-00-101xx)
30 cm, 35 cm, 40 cm
- Stiffening Bumpers:
 - 23-25 cm (KIT-00-11423)
 - 26-28 cm (KIT-00-11424)
 - 29-31 cm (KIT-00-11425)

Note: Spectra socks and heel bumpers are included in the delivery.

Complete Prosthetic Foot (incl. foot shell) Part Number

F21-00-1 **A** - -

Category: 1-9 Size: 23-30 Version: A-mid profile, S-sandal toe Side: L-left, R-right Generation: 3-regular toe, 2-sandal toe Skin Tone: L-light, M-medium, D-dark

Foot Module Part Number

F21-00-0 **A** -

Category: 1-9 Size: 23-30 Toe Configuration: RU-regular toe, Universal, SL- Sandal toe, Left, SR- Sandal toe, Right

Foot Shell Part Number

FTC- **A-1** **4-**

Generation: 2-sandal toe, 3-mid profile Size: 23-30 Skin Tone: L-light, M-medium, D-dark Toe: R-regular, S-sandal Side: L-left, R-right

See the Cosmetic Foot Shell page for more detailed ordering instructions.