

Renegade® & Renegade® LP

The industry's most innovative shock absorbing prosthetic foot

- Patented Z-Shock Technology™ is engineered to deliver a smooth gait at any speed, resulting in a comfortable walking foot that can also stand up to high impact forces
- All-carbon z-shaped heel is durable and maintenance-free
- Less than one-third the weight of alternative shock foot products, reducing socket pistoning and fatigue
- Rated to 166 kg
- 36-month warranty



Vertical Shock




Sandal Toe



EnduraCore® Technology

Renegade® & Renegade® LP

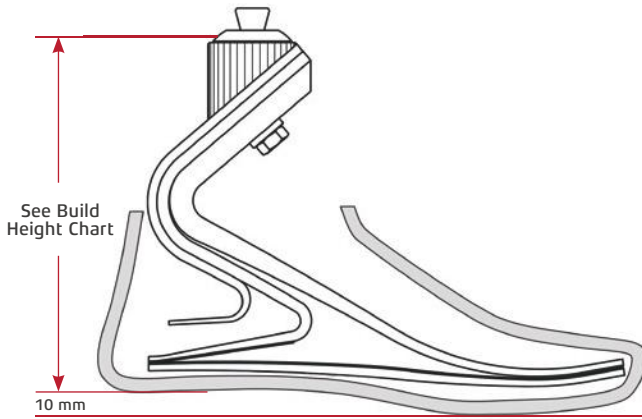

Note: Build height measurements represent a fully assembled unit.

Product ID: Renegade (RS1); Renegade LP (RS6)

Sizes: 22-31 cm

Sandal toe available for sizes 22-28 cm*

User Weight Rating: 166 kg

Heel Height: 10 mm

Connector: Male pyramid

Warranty: Foot module 36 months,
Cosmetic foot shell 6 months

*When ordering the Sandal Toe option, please contact Customer Service for lead time.

Accessories

- Pre-Cut Shin Foam: (RS1only; two piece clamshell design)
 - 22-25 cm (AFP-00-RS101)
 - 26-28 cm (AFP-00-RS102)
 - 29-31 cm (AFP-00-RS103)
- Sole Plate: (size and category specific)
 - Universal toe: (RS1-SP)
 - Sandal toe: (RS1-SP-S)
- Foot Shell Removal Tool: (ACC-00-103)
- Spectra Sock: (ACC-00-101xx)
 - 30 cm, 35 cm, 40 cm
- Stiffening Bumpers:
 - 22-25 cm (KIT-00-116)
 - 26-28 cm (KIT-00-117)
 - 29-31 cm (KIT-00-118)

Note: Spectra socks and heel bumpers are included in the delivery.

Build Height & Weight

	Size (cm)	Height (mm)	Weight (g)
Renegade	22-25	162	625
	26-28	184	765
	29-31	200	945
Renegade LP	22-25	127	625
	26-28	140	760
	29-31	149	930

Product weight includes foot shell. Measurements are approximate and will vary due to category selection and manufacturing variability.


Stiffness Category Selection

Impact Level			Weight (kg)
Low	Moderate	High	
1	1	2	44-52
1	2	3	53-59
2	3	4	60-68
3	4	5	69-77
4	5	6	78-88
5	6	7	89-100
6	7	8	101-116
7	8	9	117-130
8	9		131-147
9			148-166


Complete Prosthetic Foot (incl. foot shell) Part Number



Foot Module Part Number

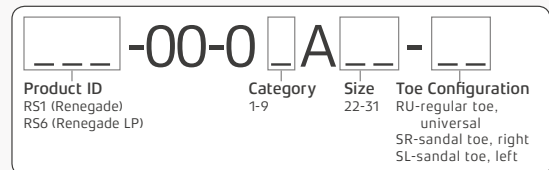

Foot Shell Part Number

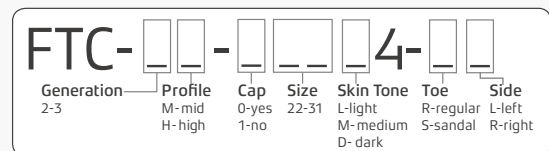

See the Cosmetic Foot Shell page for more detailed ordering instructions.