



1

2

3

4



SPECIALTY

# Silhouette<sup>®</sup> VS & Silhouette<sup>®</sup> LP-VS

## Shock absorption softens otherwise unforgiving surfaces

- Vertical shock absorption softens impact at heel strike to provide comfort and preserve residual limb health
- +/- 15 degrees of inversion/eversion provides excellent ground compliance and stability
- The lightweight, slim-profile design facilitates an easy to achieve cosmetic finish capable of pleasing the most discerning critic
- The carbon fiber sole plate with tri-durometer urethane strips facilitates a biomechanically correct ground reaction force, moving lateral to medial through the foot with each step
- The choice between 10 mm and 19 mm heel heights provides the ability to wear higher heel dress shoes or work boots while maintaining proper prosthetic alignment
- A 166 kg weight rating makes the Silhouette VS appropriate for a broad range of low to moderate impact K2 and K3 level ambulators
- 36-month warranty




Vertical Shock

Silhouette LP-VS (R15/R19)

Silhouette VS (R14/R18)

# Silhouette<sup>®</sup> VS & Silhouette<sup>®</sup> LP-VS

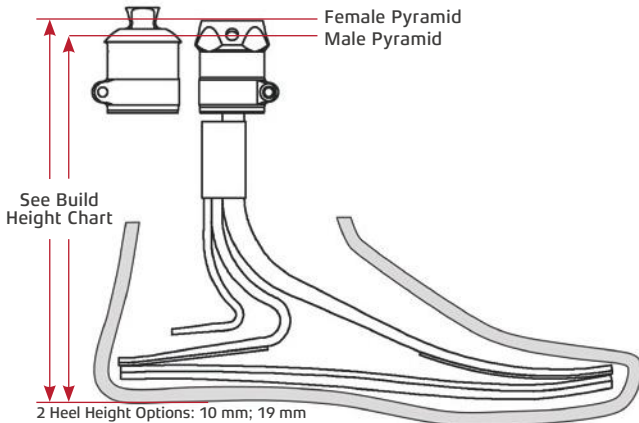

Note: Build height measurements represent a fully assembled unit.

- Product ID:** Silhouette VS (R14/R18);  
Silhouette LP-VS (R15/R19)
- Sizes:** 22-31 cm
- User Weight Rating:** 166 kg
- Heel Height:** 10 mm or 19 mm
- Connector:** Male or female pyramid, or Versa Clamp
- Warranty:** Foot module 36 months,  
Cosmetic foot shell 6 months

## Accessories

- Female pyramid, 34 mm: (ACC-00-129)
- Male pyramid, 34 mm: (ACC-00-196)
- Tube Clamp Insert: (ACC-00-124)  
(custom cut, foot module serial number required)
- Versa Clamp Adapter Kit: (KIT-00-123V) includes extended female pyramid and spacer ring
  - Extended 34 mm Female TubeClamp: (ACC-EX-123)
  - Spacer Ring: (KIT-00-123S)
- Knee Pyramid Adapter, Male: (ACC-00-195) for knees receiving a 30 mm pylons
- Replacement Sole Plate, Tri Durometer: (RS4-SP-0xAxx-xx)
- Foot Shell Removal Tool: (ACC-00-103)
- Spectra Sock: (ACC-00-101xx)  
30 cm, 35 cm, 40 cm
- Stiffening Bumpers:
  - 22-25 cm (KIT-00-12101)
  - 26-28 cm (KIT-00-12102)
  - 29-31 cm (KIT-00-12103)

Note: Spectra socks and heel bumpers are included in the delivery.

## Build Height & Weight

	Size (cm)	Height (mm)	Weight (g)
<b>Silhouette VS w/female adapter</b>	22-25	178	525
	26-28	191	645
	29-31	201	785
<b>Silhouette LP-VS w/male adapter</b>	22-25	168	430
	26-28	181	550
	29-31	191	685

Product weight includes foot shell. Measurements are approximate and will vary due to category selection and manufacturing variability.

## Stiffness Category Selection

Impact Level			Weight (kg)
Low	Moderate	High	
1	1	2	44-52
1	2	3	53-59
2	3	4	60-68
3	4	5	69-77
4	5	6	78-88
5	6	7	89-100
6	7	8	101-116
7	8	9	117-130
8	9		131-147
9			148-166

## Complete Prosthetic Foot (incl. foot shell) Part Number

   - 0    - 0    A    - M    - 3   

Product ID	Adapter	Category	Size	Side	Skin Tone
R14 (10 mm)	1-female	1-9	22-31	L-left	L-light
R15 (LP 10 mm)	2-male			R-right	M-medium
R18 (19 mm)	3-versa				D-dark
R19 (LP 19 mm)					

### Foot Module Part Number

   - 0    - 0    A    -   

Product ID	Adapter	Category	Size	Toe Configuration
R14 (10 mm)	1-female	1-9	22-31	RR-regular toe, right
R15 (LP 10 mm)	2-male			RL-regular toe, left
R18 (19 mm)	3-versa			
R19 (LP 19 mm)				

### Foot Shell Part Number

**FTC-**       -                   **4-R**   

Generation	Profile	Cap	Size	Skin Tone	Side
2-3	M-mid H-high	0-yes 1-no	22-31	L-light M-medium D-dark	L-left R-right

See the Cosmetic Foot Shell page for more detailed ordering instructions.