


Confidence & Comfort



K₁ + K₂

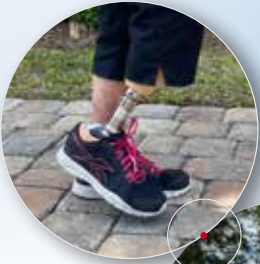

Restore™ prosthetic foot represents a breakthrough in performance for single speed ambulators (K₁ + K₂ level amputees). Featuring full length fiberglass keel in a C-shaped design, Restore provides excellent stability and support. Users will experience increased confidence when performing day-to-day activities. Prosthetists will appreciate the durable fiberglass design with easy to customize heel stiffness adjustments.

Features

- **Smooth Roll Over:** Full length fiberglass keel that provides stability and smooth roll over characteristics needed for K₁+ K₂ level ambulators.
- **Adjustable Heel Response:** A set of 3 independent heel bumpers allows the prosthetist to tune the response of the heel based on the users walking dynamics.
- **Inversion/Eversion:** Semi-split keel and the added compliance of fiberglass allow the foot to adapt to uneven terrain.
- **C-Shaped Foot Module:** C-Shape provides improved response and stability from Initial Contact through toe off, by increasing the length of the keel.

Restore™ Specifications

- Sizes: 22 cm - 30 cm
- Sandal Toe: 23 cm - 28 cm
- User Weight Rating: 136 kg
- Build Height Range: 89 mm - 99 mm
- Average Product Weight: 674 g (with footshell)
- Warranty: Foot Module 12 months (cosmetic foot shell 6 months)
- Splash and occasional water submersion
- Increased durability, flexibility and natural feel

Ordering Guide

1 Verify User's Specifications

- User Weight
- Foot Size (22 - 30 cm)
- Right or Left
- Regular Toe or Sandal Toe*


*Available Sandal Toe Sizes: 23 - 28 cm



2 Stiffness Category Selection Chart

Weight (kg)	Impact Level
45 - 59	1
60 - 81	2
82 - 109	3
110 - 136	4

3 Complete Foot Module Part Number

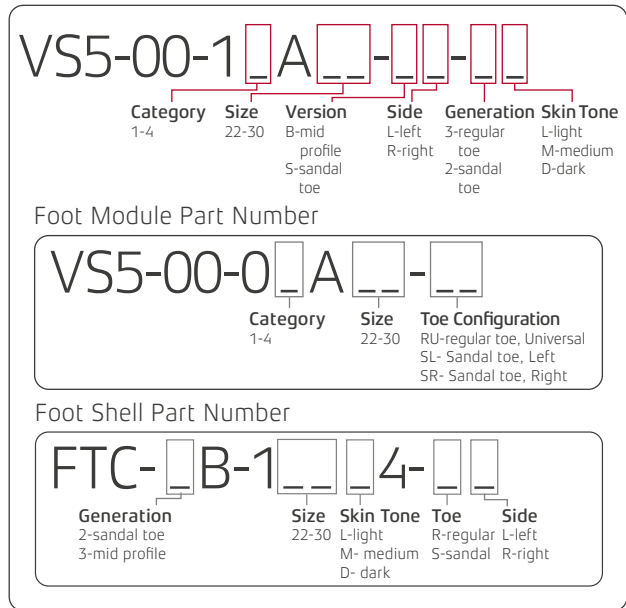

Focused Solutions

Freedom Innovations is solely focused on developing world-class lower-limb solutions in close collaboration with prosthetists and people with amputations. These solutions encompass technological innovations, service, training, educational resources and consultative support to help ensure user satisfaction.

www.freedom-innovations.eu

Freedom Innovations Europe B.V. European Head Office

Jaargetijdenweg 4
7532 SX Enschede
The Netherlands
Tel.: +31 (0) 53 20 30 300
Fax: +31 (0) 53 20 30 305

info@freedom-innovations.eu
www.freedom-innovations.eu

Global Headquarters

3 Morgan
Irvine, CA 92618
U.S.A.
Free call: +1 888 818 6777
Tel.: +1 949 672 0032
Fax: +1 949 672 0084

